

Technical specifications

Neckloop diameter: 22 cm, 9"

Weight: 46g, 1.6 oz

Cable length

BE9159: 90 cm, 3'

BE9161: 15 cm, 6"

Connectors

3.5 mm tele plug (stereo) gold plated connector, 90 degrees angle
(Break-away connector on cable)

Load impedance: 2 x 5 Ω

Magnetic output

≥ 1700 mA/m @ 15 cm, 6" distance and
2 x 50 mW input signal

Operation, storage and transport conditions

Operating temperature: 0°–35° C, 32°–95° F.

Transport and storage temperature: -10°–50° C.

Relative humidity: 15% – 90% (non condensing).

European office

Bellman & Symfon Europe AB
Södra Långebergsgatan 30,
436 32 Askim,
Sweden
Phone +46 31 68 28 20
Fax +46 31 68 28 90
E-mail info@bellman.com

Asia Pacific office

Bellman & Symfon Asia Ltd.
1401-06, Main Tower, GD Int'l Hotel
339 Huanshi Dong Rd
Guangzhou, China, 510095
Phone +86 20 8318 2950
Fax +86 20 8333 5528
E-mail info.cn@bellman.com

© 2018, Bellman & Symfon AB, all rights reserved. Bellman, Bellman & Symfon and the Bellman & Symfon logo are registered trademarks of Bellman & Symfon AB.


DESIGN FOR EARS™



bellman.com

BE9159_006MAN009

 Bellman & Symfon®
DESIGN FOR EARS

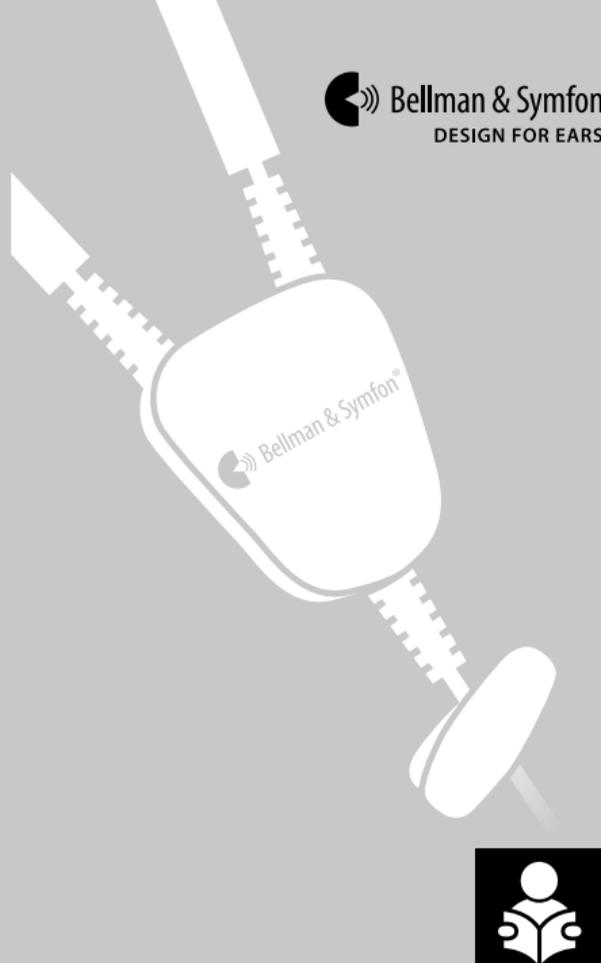


User information

- This Neckloop from Bellman & Symfon is to be used with hearing aids equipped with Telecoil.
- It is optimized for the Maxi, Sonido, Mino speech amplifiers and the Domino listening system.
- The neckloop can also be used with smart devices, TVs, radios or other equipment that features a 3.5 mm stereo headphone jack intended for low impedance headphones.
- It is discreet and comfortable, and can easily be worn under your shirt. It also incorporates a breakaway neck strap for increased safety.
- **Important!** If you have a pacemaker, we recommend you check with your general practitioner or cardiologist before using a neckloop.



Neckloop

BE9159/BE9161

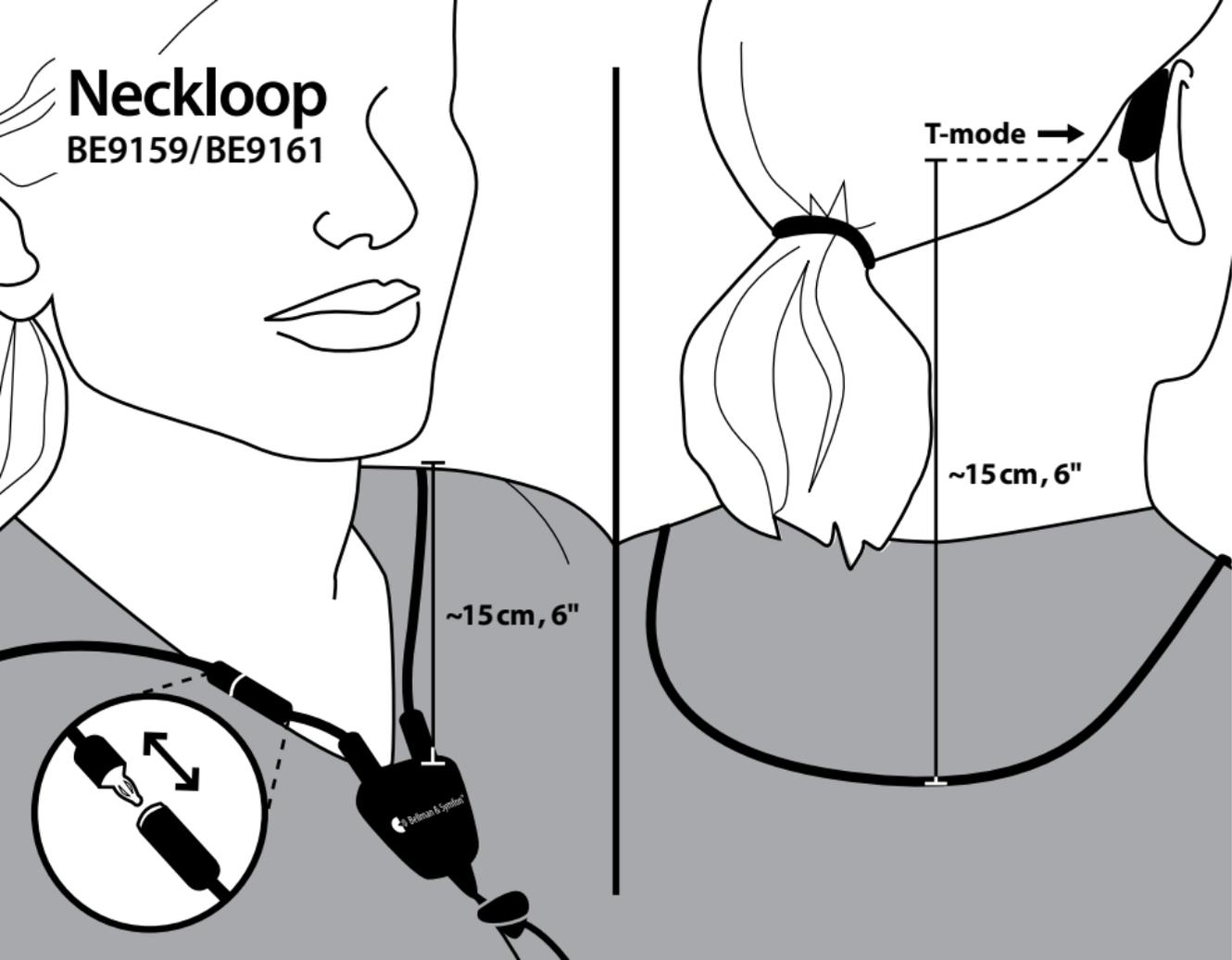

Use with Domino or Mino

